


Diode-Pumped Solid-State Blue Laser

Part No: *TECBL-30G-473*

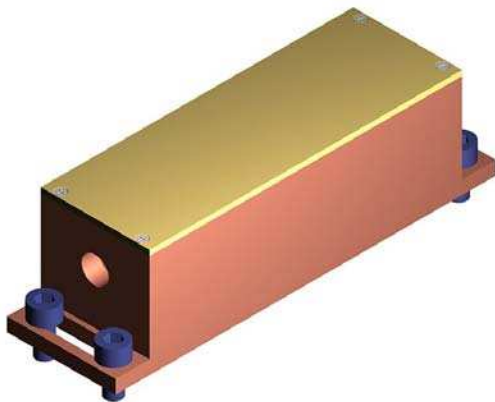

Specification

OPTICAL	
Wavelength	473 nm
Optical Output Power	30 mW
Power stability	<1%
Laser Class	Class IIIB
Laser Operation	Continuous
Laser Structure	Single Mode Laser
Spatial mode	TEM ₀₀
Beam shape	Circular (1:1.3)
Beam Divergence	<1 mrad
Spot Size at aperture at FWHM	0.3 mm
Polarization	>300:1, vertical
Warm-up time	15 min
ELECTRICAL	
Operating Voltage	5 V
Operating Current	7.5 A peak
MECHANICAL	
Laser Dimension	40 mm x 40 mm x 140 mm
Controller Dimension	69 mm x 25 mm x 153 mm
Cable Length	0.6 m
Operating Temperature	+25°C to +30°C
Storage Temperature	-40°C to +80°C
Heat Sink Requirements	5 W dissipation

Product Features

- Solid-state laser
- High Stability
- Low power noise

Applications

- Spectroscopy
- Bioanalytical
- Displays
- Light scattering

Operating Instructions

Caution note: Please read the complete instruction before attempting to operate the laser. Operating this laser using controls or performance procedures other than those specified herein may result in hazardous laser light exposure or may damage the laser.

To assemble

- Mount the laser head on the heat sink,
- Remove the shorting jumper from the laser head cable,
- Connect the laser head cable to the controller,
- Connect the power supply cable to the controller.

To start

- Plug-in the power supply to the outlet - the laser is ON,
- After 15 min warm-up time the laser is ready.

To stop the laser unplug the power supply.



Limited Warranty: One year. No warranty coverage for disassembly, modifications or damage due to abuse or misapplication.